

# STANDARD KERATOMETER




Keratometer I

*Measure in both millimeters and  
diopters, plus...unique "Split Mires"*

The Marco **KERATOMETER I** is a one position instrument for objectively measuring corneal curvature.

By simply rotating the Keratometer until a major axis is located, curvatures in both major and minor meridians can readily be measured. The Marco Keratometer assures objective accuracy in such measurements as well as in the fitting of contact lenses.

Manufactured with brass and stainless steel moving parts for durability, the Marco Keratometer incorporates such other features as precision objectives, achromatic prisms, and uniform diffuse illumination.

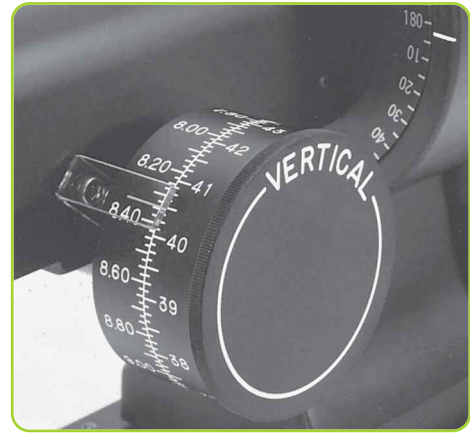

*Keratometer I with dual reading knobs measures in both mm and diopters.*

### SPECIFICATIONS & FEATURES

- Dual reading knobs measure in both millimeters and diopters.
- One position instrument. Measures both meridians without changing optical system.
- Simple vertical adjustment to fit patient.
- Two way adjustable head and chin rests.
- Dual eye-level sighting system. Facilitates horizontal alignment.
- Positive fixation. Permits rapid measurement of central corneal area.
- Coincidence focusing system. Positive measurements to 1/8-diopter accuracy.
- Wide-angle eyepiece with engraved focusing scale.
- Brass and steel moving parts assure durability.
- Precision objectives.
- Achromatic lenses.
- Fits all standard instrument trays.


**Dimensions** 16" H x 19" L x 7" W  
(406mm x 483mm x 178mm)  
Weight: 18 lb.

**Color/Finish** Black

**Standard Accessories** Dust cover, one spare bulb, 15 watt.

**Optional Accessories** Keratometer plate for instrument stand mounting.

**Stock #** 3003

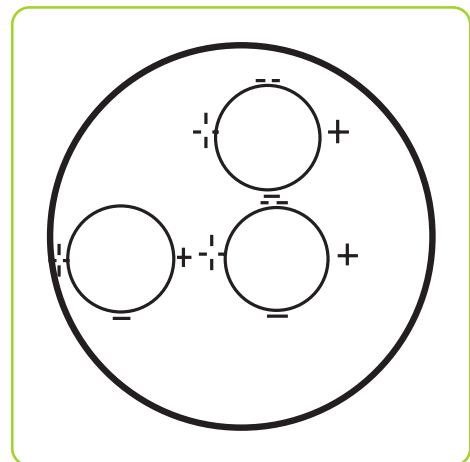

*The Keratometer I has unique "split mires" which helps to eliminate parallax.*



Marco-Certified Quality assures you of the very finest vision products, offering advanced technology, industry-leading durability and enhanced productivity for your practice. Marco-Certified Quality is your mark of confidence.



#### PLAY IT SAFE.

Protect your warranty. Buy your Marco product only from an authorized Marco distributor. Call Marco for the name of the distributor nearest you.

**T H E L E A D E R I N V I S I O N D I A G N O S T I C S ®**

